

Instructions to use the RaCAM_{IS} and RaCAM_{OD} software:

For both software, detailed instructions on how to fill each field and what they correspond to can be found in the “Instructions” tab in the menu bar of the software.

1- RaCAM_{IS}:

a. Aim

Using the RaCAM_{IS} (short for Image Segmentation) software, the user can perform image acquisition, and / or automated image processing (segmentation using ImageJ), and / or automated object recognition of the segmented images using trained CNNs. Original, segmented, and classified images can be saved. The software can also perform batch image segmentation and / or object recognition of existing images, and generate summary census data files compiling taxa counts for every sample of a core.

b. Instructions

-Open the software

-On the first row of buttons, select the operation you wish to perform (“Image acquisition” and / or “Image segmentation” and / or “Object recognition”)

Image acquisition		Image segmentation		Object recognition	
☒ Yes	☐ No	☐ Yes	☒ No	☐ Yes	☒ No

-On the second row, select the corresponding saving folder / imagej script / CNN model network_info.xml file depending on the requested operations (e.g. you don’t need to select an imagej script or CNN model if you only want to perform image acquisition; you can keep the saving folder by default).

Select saving folder	Select ImageJ script	Select CNN model
/home/esnz/Desktop	No script selected	No network_info.xml selected


-On the third row, fill in information about the core, sample and magnification you are using. These can account for other metadata if you are not using a core and sample format. If the magnification is not in the dropdown list, you can directly write it (providing magnification has no impact on any data processing, it is just a metadata that will be written in the filename).

Core name	Sample name	Magnification
Core_name	Sample_name	x10

-On the fourth, fifth and sixth rows, adjust the parameters as you need. More information about each parameter can be found in the “Instructions” tab directly on the software.

Brightness		Contrast	Exposure time	
0.05		1.3	500	
Image rotation		Original image quality	Image type	
180		100	☒ Greyscale ☐ Colour	
Image object quality		Min. object size filter	CNN threshold	
90		Script default	0.3	

-The “Save profile” and “Load profile” can be used to save all these parameters before shutting down the software. Only one profile can be saved.

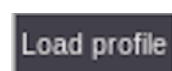

-Clicking on the “Live preview” button shows a livestream preview from your microscope so you can adjust the parameters, focal depth and slide position before capturing and saving an image. The preview window needs to be closed before using any other function.



-Clicking on the “Snapshot” button captures an image and saves it in the “saving folder” (again, we recommend keeping the saving folder by default on the Desktop). The image will be saved in the “/Desktop/RaCAM_IS_output/Image_acquisition/” directory following the core name and sample name provided. If you also selected “Image segmentation”: “Yes”, and / or “Image recognition”: “Yes”, the corresponding images will also be saved in the “/Desktop/RaCAM_IS_output/Image_processing/” and “/Desktop/RaCAM_IS_output /Image_recognition/” directories, respectively.



-Clicking on the “Batch processing” button allows you to process an existing batch of image (either do segmentation on a batch of existing original FOV images, or identify a batch of segmented images, or both) by unselecting “Image_acquisition” and selecting “Image_processing” and / or “Image_recognition”: “Yes”. You will be asked to choose and enter a folder, which should be a “core” directory containing “sample” directories, containing images.



-Clicking the “Census data” button generates a census table containing taxa in columns and samples in rows, by counting the number of images located in “class” subdirectories in “sample” directories in the “/Desktop/RaCAM_IS_output/Temp/” directory automatically generated when automated recognition is performed. After performing the census counting, this directory is automatically deleted.

Census data

2- RaCAM_{OD}:

a. Aim

Using the RaCAM_{OD} (short for Object Detection) software, the user can perform image acquisition, and automated object detection directly on captured field of views using a trained .onnx YOLO model. Original and annotated field of view images can be saved, as well as census data txt files for each sample. The software can also perform batch object detection for existing field of view images, and generate summary census data files compiling taxa counts for every sample of a core. The user can also perform live stream object detection using the same .onnx YOLO model. The software is also capable of automatically capture images when a target species is detected.

b. Instructions

-Open the software

-On the first row of buttons, select the operation you wish to perform: “Image acquisition” and / or “Object detection” and / or “Show annotated images” (that will automatically open the annotated images after automated acquisition and object detection for manual control).

Image acquisition		Object detection		Show annotated images	
☒ Yes	☐ No	☐ Yes	☒ No	☐ Yes	☒ No

-On the second row, select the corresponding saving folder and / or .onnx model file and confidence threshold for object detection depending on the requested operations (e.g. you don’t need to select an .onnx model if you only want to perform image acquisition; you can keep the saving folder by default).

Select saving folder	Select .onnx model	Confidence threshold
/home/esnz/Desktop	No .onnx model selected	0.3

-On the third row, fill in informations about the core, sample and magnification you are using. These can account for other metadata if you are not using a core and sample format. If the magnification is not in the dropdown list, you can directly write it (providing magnification has no impact on any data processing, it is just a metadata that will be written in the filename).

Core name	Sample name	Magnification
Core_name	Sample_name	x10

-On the fourth and fifth rows, adjust the parameters as you need. More information about each parameter can be found in the “Instructions” tab directly on the software.

Brightness	Contrast	Exposure time
0.05	1.3	500
Image rotation	Original image quality	Image type
180	100	☒ Greyscale ☐ Colour

-The “Save profile” and “Load profile” can be used to save all these parameters before shutting down the software. Only one profile can be saved.

Load profile

-Clicking on the “Live preview” button shows a livestream preview from your microscope so you can adjust the parameters, focal depth and slide position before capturing and saving an image. The preview window needs to be closed before using any other function.

Live preview

-By selecting “Image acquisition”: “No” and “Object detection”: “Yes”, the “Live preview” button will change to “Live detection”, allowing you to perform automated live detection using a selected .onnx model. A pop-up window will ask you for a target species: if you write a target class name exactly as learnt by the model, a snapshot will automatically be taken by the software every time this target class is detected with a cooldown of 10 seconds, allowing to manually drag the microscopic slide under the microscope and let the system capture images when needed. If you leave this field empty, no image will be captured automatically.

-Clicking on the “Snapshot” button captures an image and saves it in the “saving folder” (again, we recommend keeping the saving folder by default on the Desktop). The image will be saved in the “/Desktop/RaCAM_OD_output/Image_acquisition/” directory following the core name and sample name provided. If you also select “Object

detection”: “Yes”, the corresponding images will be saved in the “/Desktop/RaCAM_OD_output/Image_recognition/” directory.

Snapshot

-Clicking on the “Batch processing” button allow you to process an existing batch of image (do object detection on a batch of existing original FOV images) by unselecting “Image_acquisition” and selecting “Object detection”: “Yes. You will be asked to choose and enter a folder, which should be a “core” directory containing “sample” directories, containing images.

Batch processing

-Clicking the “Census data” button will generate a census table containing taxa in columns and samples in rows, by compiling the individual .txt files generated during the object detection step for each FOV image and located in the “/Desktop/RaCAM_OD_output/Raw_census_data/” directory.

Census data